


Handlin

Manipolatore/iniettore uterino
Uterine manipulator/Injector

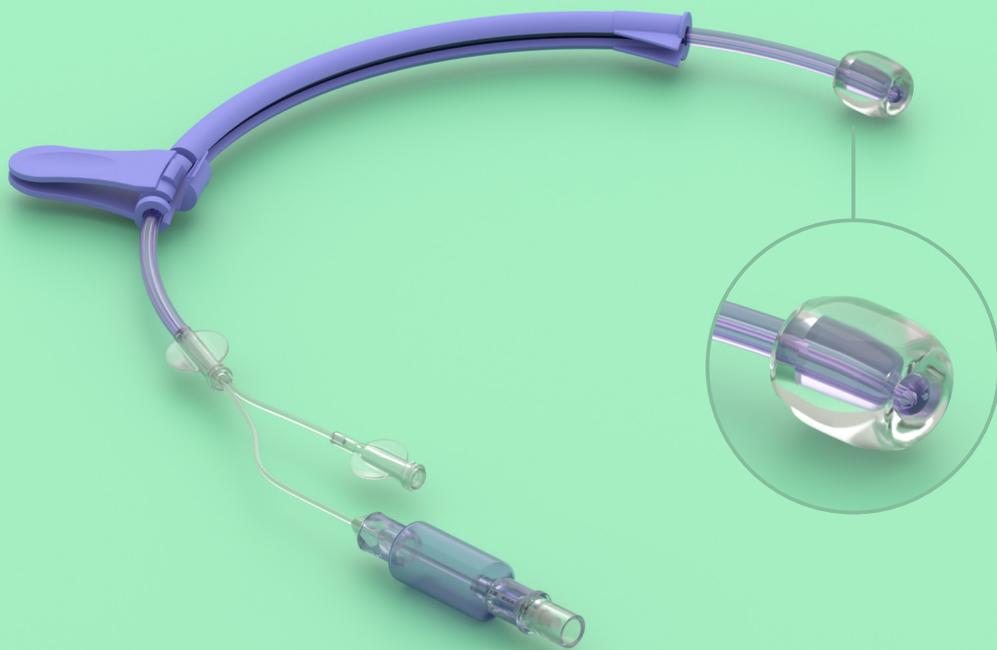

Impiegato in procedure dove viene richiesta la manipolazione dell'utero e/o per iniettare liquidi di contrasto (quali laparoscopia diagnostica, esami di fertilità, etc.).

È caratterizzato da una forma ricurva che ne facilita l'inserimento.

Provvisto di:

- un catetere interno a due lumi con scala di misura da 4 a 10 cm (risoluzione 1 cm)* per determinare la corretta profondità di inserimento;
- un'impugnatura rigida;
- un sistema di bloccaggio semplice e sicuro che impedisce movimenti accidentali durante la procedura;
- un palloncino in PVC da 10 cc e 20 mm di diametro alloggiato sulla punta del catetere in modo da minimizzare la possibilità di perforazione accidentale ed evitare il reflusso del mezzo di contrasto;
- una valvola unidirezionale per gonfiare il palloncino, assemblata a un polmoncino di riferimento che permette di visualizzare la reale condizione del palloncino interno;
- un raccordo luer lock per l'iniezione del mezzo di contrasto.

* Fornisce una misura approssimativa.




It is used for procedures where uterus manipulation and/or injection of contrast media is requested (e.g. diagnostic laparoscopy, fertility examinations, etc.). It has a curved shape for an easy insertion.

It is composed by:

- double lumen inner catheter with graduation from 4 to 10 cm (1 cm resolution)* to assess correct insertion depth;
- stiff handle allowing a firm and safe manipulation during the procedure;
- simple locking system to avoid accidental movements during the procedure;
- 10 cc PVC balloon with diameter 20 mm placed on the catheter tip to minimize any accidental perforation and avoid contrast media reflux;
- one-way inflation valve directly connected to a pilot balloon showing the real condition of the inner balloon;
- luer lock connector for contrast media injection.

* It gives an approximate measurement.

Cod.	Descrizione / Description	Pz/Scatola - Pcs/Box
720199	Ø esterno 5 mm O.D 5 mm	10 pz./pcs.

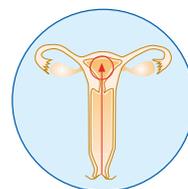

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked