


Rigidinjector

Iniettore uterino per isterosalpingografia
Hysterosalpingography uterine injector

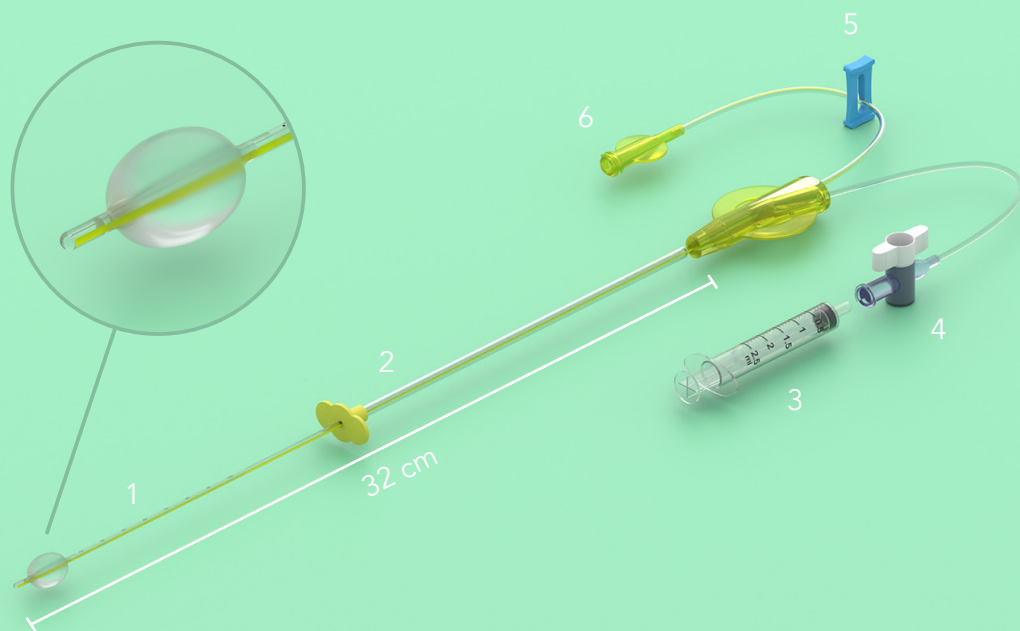

 Per l'iniezione di mezzi di contrasto (radiopaco / blu di metilene) nelle procedure di isterosalpingografia. Sono costituiti da:


- un introduttore rigido dotato di fermo scorrevole (2);
- un catetere rigido radiopaco a 2 lumi (1) dotato di un palloncino ovale in gomma termoplastica da 2,5 cc sulla parte distale che assicura un alloggiamento stabile del catetere;
- rubinetto (4) per regolare il gonfiaggio del palloncino;
- clamp stringitubo per evitare il reflusso di liquido (5);
- un raccordo luer lock femmina (6) per l'iniezione del mezzo di contrasto;
- una siringa da 2,5 cc (3).

 To inject contrast media (radio-opaque / methylene blue) in hysterosalpingography diagnostic procedure.

The catheters are composed by:

- rigid introducer provided with a mobile positioner (2);
- double lumen radiopaque rigid catheter (1) with a thermoplastic rubber oval balloon of 2.5 cc on the distal tip, that ensures a firm positioning of the catheter;
- stopcock (4) to regulate balloon inflation;
- tubing clamp to avoid the liquid backflow (5);
- female luer-lock connector (6) for contrast media injection;
- a 2.5 cc syringe (3).

Cod.	Descrizione / Description	Pz/Scatola - Pcs/Box
720239	Ø esterno 1,67 mm (5 Fr) O.D 1.67 mm (5 Fr)	10 pz./pcs.
720225	Ø esterno 2 mm (6 Fr) O.D 2 mm (6 Fr)	
720251	Ø esterno 2,33 mm (7 Fr) O.D 2.33 mm (7 Fr)	

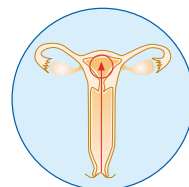

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked